

FUCHSIN BASIC

Safety Data Sheet

according to Regulation (EU) 2015/830

Date of issue: 16/02/2017 Version: 0.0


Doc. No: SDS-924.D02/1



SECTION 1: Identification of the substance/mixture and of the company/undertaking

1.1. Product identifier

Product form : Substance
Trade name : FUCHSIN BASIC FOR MICROSCOPY
EC No : 211-189-6
CAS No : 632-99-5
Type of product : Pure substance
Formula : C₂₀H₁₉N₃.HCl
Chemical structure :

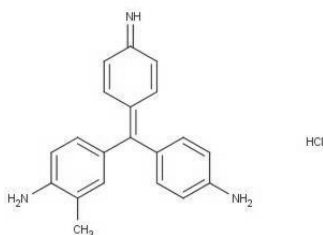

Synonyms

: (4-(4-aminophenyl)(4-iminocyclohexa-2,5-dienylidene)methyl)-2-methylaniline hydrochloride / (4-amino-3-methylphenyl)bis(4-aminophenyl)methanol / 12418 Red / 4-((4-aminophenyl)(4-imino-2,5-cyclohexadien-1-ylidene)methyl)-2-methylbenzenamine HCl / Aizen magenta / aniline red / astra fuchsine B / basic fuchsin (=C.I. basic violet 14) / basic fuchsine (=C.I. basic violet 14) / basic magenta (=C.I. basic violet 14) / basic magenta E-200 (=C.I. basic violet 14) / basic violet 14 / benzenamine, 4-((4-aminophenyl)(4-imino-2,5-cyclohexadien-1-ylidene)methyl)-2-methyl-, monohydrochloride / C.I. 42510 / C.I. basic violet 14, monohydrochloride / C.I. solvent red 41 / calcozine fuchsine HO / calcozine magenta RIN / calcozine magenta RTN / calcozine magenta XX / cerise B / C-WR violet 8 / diabasic magenta / diamond fuchsine (=C.I. basic violet 14) / fuchsin (=C.I. basic violet 14) / fuchsine (=C.I. basic violet 14) / fuchsine A / fuchsine CS / fuchsine G / fuchsine HO / fuchsine N / fuchsine RTN / fuchsine SBP / fuchsine Y / magenta DP / magenta E / magenta G / magenta I / magenta PN / magenta powder N / magenta S / magenta Super Fine / methyl fuchsine / orient basic magenta / rosaniline chloride / rosaniline hydrochloride / rosanilinium chloride / rosanilinium hydrochloride / violet 14, basic / violet zasadita 14

1.2. Relevant identified uses of the substance or mixture and uses advised against

1.2.1. Relevant identified uses

Industrial/Professional use spec : Industrial
For professional use only
Use of the substance/mixture : Dyestuff/pigment
Laboratory chemical

1.2.2. Uses advised against

No additional information available

1.3. Details of the supplier of the safety data sheet

ISOLAB Laborgeräte GmbH
Am Dillhof 2 - 63863 Eschau / GERMANY
Tel: +49 93 74 / 978 55-0
Fax: +49 93 74 / 978 55-29
prodsafe@isolab.de

1.4. Emergency telephone number

Country	Organisation/Company	Address	Emergency number	Comment
Germany	Giftnotruf der Charité CBF, Haus VIII (Wirtschaftsgebäude), UG	Hindenburgdamm 30 12203 Berlin	+49 30 19240	

SECTION 2: Hazards identification

2.1. Classification of the substance or mixture

Classification according to Regulation (EC) No. 1272/2008 [CLP]

Acute Tox. 4 (Oral) H302

STOT SE 2 H371

FUCHSIN BASIC

Safety Data Sheet

according to Regulation (EU) 2015/830

Date of issue: 2/16/2017

Version: 0.0

Doc. No: SDS-924.D02/1



Full text of hazard classes and H-statements : see section 16

Adverse physicochemical, human health and environmental effects

Suspected of causing cancer. Harmful if swallowed. May cause respiratory irritation. Causes skin irritation. Causes serious eye irritation.

2.2. Label elements

Labelling according to Regulation (EC) No. 1272/2008 [CLP]

Hazard pictograms (CLP) :



GHS07

GHS08

Signal word (CLP) :

Warning

Hazard statements (CLP) :

H302 - Harmful if swallowed
H371 - May cause damage to organs.

Precautionary statements (CLP) :

P260 - Do not breathe dust, fume, gas, mist, vapours, spray.
P264 - Wash thoroughly after handling.
P270 - Do not eat, drink or smoke when using this product.
P301+P312 - IF SWALLOWED : Call a POISON CENTER or doctor if you feel unwell.
P309+P311 - IF exposed or if you feel unwell : Call a POISON CENTER or doctor.
P405 - Store locked up.
P501 - Dispose of this material and its container to hazardous or special waste collection point, in accordance with local, regional, national and/or international regulation.

2.3. Other hazards

This substance/mixture does not meet the PBT criteria of REACH regulation, annex XIII

This substance/mixture does not meet the vPvB criteria of REACH regulation, annex XIII

SECTION 3: Composition/information on ingredients

3.1. Substances

Name	Product identifier	%
FUCHSIN BASIC	(CAS No) 632-99-5 (EC No) 211-189-6	100

Full text of H-statements: see section 16

3.2. Mixtures

Not applicable

SECTION 4: First aid measures

4.1. Description of first aid measures

First-aid measures general	: Check the vital functions. Unconscious: maintain adequate airway and respiration. Respiratory arrest: artificial respiration or oxygen. Cardiac arrest: perform resuscitation. Victim conscious with laboured breathing: half-seated. Victim in shock: on his back with legs slightly raised. Vomiting: prevent asphyxia/aspiration pneumonia. Prevent cooling by covering the victim (no warming up). Keep watching the victim. Give psychological aid. Keep the victim calm, avoid physical strain. Depending on the victim's condition: doctor/hospital. IF exposed or concerned: Get medical advice/attention. Call a poison center or a doctor if you feel unwell.
First-aid measures after inhalation	: Remove person to fresh air and keep comfortable for breathing. Remove the victim into fresh air. Respiratory problems: consult a doctor/medical service. Call a poison center or a doctor if you feel unwell.
First-aid measures after skin contact	: Wash immediately with lots of water. Soap may be used. Take victim to a doctor if irritation persists. Wash skin with plenty of water. Take off contaminated clothing. If skin irritation occurs: Get medical advice/attention.
First-aid measures after eye contact	: Rinse immediately with plenty of water. Take victim to an ophthalmologist. Rinse cautiously with water for several minutes. Remove contact lenses, if present and easy to do. Continue rinsing. If eye irritation persists: Get medical advice/attention.
First-aid measures after ingestion	: Rinse mouth with water. Immediately after ingestion: give lots of water to drink. Do not induce vomiting. Call Poison Information Centre (www.big.be/antigif.htm). Consult a doctor/medical service if you feel unwell. Ingestion of large quantities: immediately to hospital. Doctor: administration of chemical antidote. Rinse mouth. Call a poison center or a doctor if you feel unwell.

FUCHSIN BASIC

Safety Data Sheet

according to Regulation (EU) 2015/830

Date of issue: 2/16/2017

Version: 0.0

Doc. No: SDS-924.D02/1



4.2. Most important symptoms and effects, both acute and delayed

Symptoms/injuries after inhalation	: AFTER INHALATION OF DUST: Dry/sore throat. May cause respiratory irritation.
Symptoms/injuries after skin contact	: Red skin. Tingling/irritation of the skin. Irritation.
Symptoms/injuries after eye contact	: Irritation of the eye tissue. Lacrimation. Redness of the eye tissue. Eye irritation.
Symptoms/injuries after ingestion	: Abdominal pain. Nausea. Vomiting.
Chronic symptoms	: ON CONTINUOUS/REPEATED EXPOSURE/CONTACT: Skin rash/inflammation.

4.3. Indication of any immediate medical attention and special treatment needed

Treat symptomatically.

SECTION 5: Fire-fighting measures

5.1. Extinguishing media

Suitable extinguishing media	: Water spray. Polyvalent foam. Polymer foam. ABC powder. Carbon dioxide. Sand/earth. Dry powder.
Unsuitable extinguishing media	: No unsuitable extinguishing media known.

5.2. Special hazards arising from the substance or mixture

Fire hazard	: DIRECT FIRE HAZARD. No data available on direct fire hazard. INDIRECT FIRE HAZARD. No data available on indirect fire hazard.
Explosion hazard	: DIRECT EXPLOSION HAZARD. Fine dust is explosive with air. INDIRECT EXPLOSION HAZARD. Dust cloud can be ignited by a spark.
Hazardous decomposition products in case of fire	: Toxic fumes may be released.

5.3. Advice for firefighters

Precautionary measures fire	: Exposure to fire/heat: keep upwind. Exposure to fire/heat: consider evacuation. Exposure to fire/heat: have neighbourhood close doors and windows.
Firefighting instructions	: Cool tanks/drums with water spray/remove them into safety. Dilute toxic gases with water spray. Take account of environmentally hazardous firefighting water. Use water moderately and if possible collect or contain it.
Protection during firefighting	: Heat/fire exposure: compressed air/oxygen apparatus. Do not attempt to take action without suitable protective equipment. Self-contained breathing apparatus. Complete protective clothing.

SECTION 6: Accidental release measures

6.1. Personal precautions, protective equipment and emergency procedures

6.1.1. For non-emergency personnel

Protective equipment	: Gloves. Face-shield. Protective clothing. Dust cloud production: compressed air/oxygen apparatus. Dust cloud production: dust-tight suit.
Emergency procedures	: Ventilate spillage area. Mark the danger area. Prevent dust cloud formation. No naked flames. Wash contaminated clothes. Avoid breathing dust/fume/gas/mist/vapours/spray. Avoid contact with skin and eyes.
Measures in case of dust release	: In case of dust production: keep upwind. Dust production: have neighbourhood close doors and windows. Dust production: stop engines and no smoking. In case of dust production: no naked flames or sparks. Dust: spark-/explosionproof appliances/lighting equipment.

6.1.2. For emergency responders

Protective equipment	: Do not attempt to take action without suitable protective equipment. For further information refer to section 8: "Exposure controls/personal protection".
----------------------	---

6.2. Environmental precautions

Avoid release to the environment. Prevent soil and water pollution. Prevent spreading in sewers.

6.3. Methods and material for containment and cleaning up

For containment	: Contain released substance, pump into suitable containers. Plug the leak, cut off the supply. Dam up the solid spill. Knock down/dilute dust cloud with water spray. Powdered form: no compressed air for pumping over spills.
Methods for cleaning up	: Mechanically recover the product. Prevent dust cloud formation. Carefully collect the spill/leftovers. Collect the spill only if it is in a dry state. Scoop solid spill into closing containers. Powdered: do not use compressed air for pumping over spills. Clean contaminated surfaces with an excess of water. Take collected spill to manufacturer/competent authority. Wash clothing and equipment after handling. Notify authorities if product enters sewers or public waters.
Other information	: Dispose of materials or solid residues at an authorized site.

FUCHSIN BASIC

Safety Data Sheet

according to Regulation (EU) 2015/830

Date of issue: 2/16/2017

Version: 0.0

Doc. No: SDS-924.D02/1



6.4. Reference to other sections

For further information refer to section 13.

SECTION 7: Handling and storage

7.1. Precautions for safe handling

Precautions for safe handling

: Comply with the legal requirements. Remove contaminated clothing immediately. Clean contaminated clothing. Do not discharge the waste into the drain. Powdered form: no compressed air for pumping over. Avoid raising dust. Keep away from naked flames/heat. Finely divided: spark- and explosionproof appliances. Finely divided: keep away from ignition sources/sparks. Observe strict hygiene. Keep container tightly closed. Carry operations in the open/under local exhaust/ventilation or with respiratory protection. Obtain special instructions before use. Do not handle until all safety precautions have been read and understood. Wear personal protective equipment. Use only outdoors or in a well-ventilated area. Avoid breathing dust/fume/gas/mist/vapours/spray. Avoid contact with skin and eyes.

Hygiene measures

: Wash contaminated clothing before reuse. Do not eat, drink or smoke when using this product. Always wash hands after handling the product.

7.2. Conditions for safe storage, including any incompatibilities

Storage conditions

: Store locked up. Store in a well-ventilated place. Keep container tightly closed. Keep cool.

Heat and ignition sources

: KEEP SUBSTANCE AWAY FROM: heat sources. ignition sources.

Information on mixed storage

: KEEP SUBSTANCE AWAY FROM: oxidizing agents. (strong) acids. (strong) bases. water/moisture.

Storage area

: Store in a cool area. Store in a dry area. Keep container in a well-ventilated place. Meet the legal requirements.

Special rules on packaging

: SPECIAL REQUIREMENTS: closing. correctly labelled. meet the legal requirements. Secure fragile packagings in solid containers.

Packaging materials

: SUITABLE MATERIAL: No data available. MATERIAL TO AVOID: No data available.

7.3. Specific end use(s)

No additional information available

SECTION 8: Exposure controls/personal protection

8.1. Control parameters

No additional information available

8.2. Exposure controls

Appropriate engineering controls:

Ensure good ventilation of the work station.

Personal protective equipment:

Gloves.

Materials for protective clothing:

GIVE EXCELLENT RESISTANCE: No data available. GIVE GOOD RESISTANCE: No data available. GIVE LESS RESISTANCE: No data available. GIVE POOR RESISTANCE: No data available

Hand protection:

Gloves

Eye protection:

Face shield. In case of dust production: protective goggles. Safety glasses

Type	Use	Characteristics	Standard
Face shield, Safety glasses	Dust		

Skin and body protection:

Protective clothing. In case of dust production: head/neck protection. In case of dust production: dustproof clothing

Respiratory protection:

Dust production: dust mask with filter type P2

Device	Filter type	Condition	Standard
Dust mask	Type P2	Dust protection	

FUCHSIN BASIC

Safety Data Sheet

according to Regulation (EU) 2015/830

Date of issue: 2/16/2017

Version: 0.0

Doc. No: SDS-924.D02/1



Environmental exposure controls:

Avoid release to the environment.

SECTION 9: Physical and chemical properties

9.1. Information on basic physical and chemical properties

Physical state	: Solid
Appearance	: Crystalline solid.
Molecular mass	: 337.85 g/mol
Colour	: Dark green.
Odour	: Odourless.
Odour threshold	: No data available
pH	: No data available
Relative evaporation rate (butylacetate=1)	: No data available
Melting point	: 268-270 °C
Freezing point	: Not applicable
Boiling point	: No data available
Flash point	: Not applicable
Auto-ignition temperature	: Not applicable
Decomposition temperature	: 200 °C
Flammability (solid, gas)	: Non flammable
Vapour pressure	: No data available
Relative vapour density at 20 °C	: No data available
Relative density	: Not applicable
Density	: 1220 kg/m ³
Solubility	: No data available
Log Pow	: 2.92 (Estimated value)
Viscosity, kinematic	: Not applicable
Viscosity, dynamic	: No data available
Explosive properties	: No data available
Oxidising properties	: No data available
Explosive limits	: Not applicable

9.2. Other information

VOC content	: 0 %
Other properties	: Substance has acid reaction.

SECTION 10: Stability and reactivity

10.1. Reactivity

On burning: release of toxic and corrosive gases/vapours (nitrous vapours, hydrogen chloride, carbon monoxide - carbon dioxide).

10.2. Chemical stability

Stable under normal conditions.

10.3. Possibility of hazardous reactions

No dangerous reactions known under normal conditions of use.

10.4. Conditions to avoid

None under recommended storage and handling conditions (see section 7).

10.5. Incompatible materials

No additional information available

10.6. Hazardous decomposition products

Under normal conditions of storage and use, hazardous decomposition products should not be produced.

FUCHSIN BASIC

Safety Data Sheet

according to Regulation (EU) 2015/830

Date of issue: 2/16/2017

Version: 0.0

Doc. No: SDS-924.D02/1



SECTION 11: Toxicological information

11.1. Information on toxicological effects

Acute toxicity	: Oral: Harmful if swallowed.
Skin corrosion/irritation	: Causes skin irritation. pH: 5 - 6 (0.1 %)
Serious eye damage/irritation	: Causes serious eye irritation. pH: 5 - 6 (0.1 %)
Respiratory or skin sensitisation	: Not classified
Germ cell mutagenicity	: Not classified
Carcinogenicity	: Suspected of causing cancer.
Reproductive toxicity	: Not classified
STOT-single exposure	: May cause respiratory irritation.
STOT-repeated exposure	: Not classified
Aspiration hazard	: Not classified
IARC group	: 2B

SECTION 12: Ecological information

12.1. Toxicity

Ecology - general	: The product is not considered harmful to aquatic organisms nor to cause long-term adverse effects in the environment.
Ecology - air	: Not classified as dangerous for the ozone layer (Regulation (EC) No 1005/2009).
Ecology - water	: Severe water pollutant (surface water). No data available on ecotoxicity.

12.2. Persistence and degradability

FUCHSIN BASIC (632-99-5)	
Persistence and degradability	Biodegradability in water: no data available.

12.3. Bioaccumulative potential

FUCHSIN BASIC (632-99-5)	
Log Pow	2.92 (Estimated value)
Bioaccumulative potential	Low potential for bioaccumulation (Log Kow < 4).

12.4. Mobility in soil

No additional information available

12.5. Results of PBT and vPvB assessment

FUCHSIN BASIC (632-99-5)	
This substance/mixture does not meet the PBT criteria of REACH regulation, annex XIII	
This substance/mixture does not meet the vPvB criteria of REACH regulation, annex XIII	

12.6. Other adverse effects

No additional information available

SECTION 13: Disposal considerations

13.1. Waste treatment methods

Waste treatment methods	: Dispose of contents/container in accordance with licensed collector's sorting instructions.
Product/Packaging disposal recommendations	: Remove to an authorized plant for the destruction, neutralization and elimination of hazardous waste. Remove waste in accordance with local and/or national regulations. Hazardous waste shall not be mixed together with other waste. Different types of hazardous waste shall not be mixed together if this may entail a risk of pollution or create problems for the further management of the waste. Hazardous waste shall be managed responsibly. All entities that store, transport or handle hazardous waste shall take the necessary measures to prevent risks of pollution or damage to people or animals.
Additional information	: LWCA (the Netherlands): KGA category 06. Can be considered as non hazardous waste according to Directive 2008/98/EC.
European List of Waste (LoW) code	: 07 03 99 - wastes not otherwise specified

FUCHSIN BASIC

Safety Data Sheet

according to Regulation (EU) 2015/830

Date of issue: 2/16/2017

Version: 0.0

Doc. No: SDS-924.D02/1



SECTION 14: Transport information

In accordance with ADR / RID / IMDG / IATA / ADN

ADR	IMDG	IATA	ADN	RID
14.1. UN number				
Not applicable	Not applicable	Not applicable	Not applicable	Not applicable
14.2. UN proper shipping name				
Not applicable	Not applicable	Not applicable	Not applicable	Not applicable
Not applicable	Not applicable	Not applicable	Not applicable	Not applicable
14.3. Transport hazard class(es)				
Not applicable	Not applicable	Not applicable	Not applicable	Not applicable
Not applicable	Not applicable	Not applicable	Not applicable	Not applicable
14.4. Packing group				
Not applicable	Not applicable	Not applicable	Not applicable	Not applicable
14.5. Environmental hazards				
Not applicable	Not applicable	Not applicable	Not applicable	Not applicable
No supplementary information available				

14.6. Special precautions for user

- Overland transport

Not applicable

- Transport by sea

Not applicable

- Air transport

Not applicable

- Inland waterway transport

Not applicable

- Rail transport

Not applicable

14.7. Transport in bulk according to Annex II of Marpol and the IBC Code

Not applicable

SECTION 15: Regulatory information

15.1. Safety, health and environmental regulations/legislation specific for the substance or mixture

15.1.1. EU-Regulations

No REACH Annex XVII restrictions

FUCHSIN BASIC is not on the REACH Candidate List

FUCHSIN BASIC is not on the REACH Annex XIV List

VOC content : 0 %

15.1.2. National regulations

Germany

VwVwS Annex reference : Water hazard class (WGK) 3, severe hazard to waters (Classification according to VwVwS, Annex 1 or 2; ID No. 857)

WGK remark : Classification water polluting in compliance with Verwaltungsvorschrift wassergefährdender Stoffe (VwVwS) of 27 July 2005 (Anhang 2)

12th Ordinance Implementing the Federal Immission Control Act - 12.BImSchV : Is not subject of the 12. BImSchV (Hazardous Incident Ordinance)

Netherlands

SZW-lijst van kankerverwekkende stoffen : The substance is not listed

SZW-lijst van mutagene stoffen : The substance is not listed

FUCHSIN BASIC

Safety Data Sheet

according to Regulation (EU) 2015/830

Date of issue: 2/16/2017

Version: 0.0

Doc. No: SDS-924.D02/1



NIET-limitatieve lijst van voor de voortplanting giftige stoffen – Borstvoeding : The substance is not listed

NIET-limitatieve lijst van voor de voortplanting giftige stoffen – Vruchtbaarheid : The substance is not listed

NIET-limitatieve lijst van voor de voortplanting giftige stoffen – Ontwikkeling : The substance is not listed

Denmark

Recommendations Danish Regulation : Young people below the age of 18 years are not allowed to use the product
Pregnant/breastfeeding women working with the product must not be in direct contact with the product
The requirements from the Danish Working Environment Authorities regarding work with carcinogens must be followed during use and disposal

15.2. Chemical safety assessment

No chemical safety assessment has been carried out

SECTION 16: Other information

Abbreviations and acronyms:

ADN	European Agreement concerning the International Carriage of Dangerous Goods by Inland Waterways
ADR	European Agreement concerning the International Carriage of Dangerous Goods by Road
ATE	Acute Toxicity Estimate
BCF	Bioconcentration factor
CLP	Classification Labelling Packaging Regulation; Regulation (EC) No 1272/2008
EC50	Median effective concentration
IATA	International Air Transport Association
IMDG	International Maritime Dangerous Goods
LC50	Median lethal concentration
LD50	Median lethal dose
PBT	Persistent Bioaccumulative Toxic
REACH	Registration, Evaluation, Authorisation and Restriction of Chemicals Regulation (EC) No 1907/2006
RID	Regulations concerning the International Carriage of Dangerous Goods by Rail
SDS	Safety Data Sheet
vPvB	Very Persistent and Very Bioaccumulative

Data sources : REGULATION (EC) No 1272/2008 OF THE EUROPEAN PARLIAMENT AND OF THE COUNCIL of 16 December 2008 on classification, labelling and packaging of substances and mixtures, amending and repealing Directives 67/548/EEC and 1999/45/EC, and amending Regulation (EC) No 1907/2006.

Full text of H- and EUH-statements:

Acute Tox. 4 (Oral)	Acute toxicity (oral), Category 4
STOT SE 2	Specific Target Organ Toxicity - Single exposure - Category 2
H302	Harmful if swallowed
H371	May cause damage to organs.

SDS ISOLAB

This information is based on our current knowledge and is intended to describe the product for the purposes of health, safety and environmental requirements only. It should not therefore be construed as guaranteeing any specific property of the product